



## Contents

- 1. Comparison: Exterior
- 2. Platform Design
- 3. Wall Thickness
- 4. Body Design
- 5. Surgical Sequence
- 6. Apex
- 7. Result



## Comparison: Exterior





A BioSeal & Microgroove
Wall Thickness

Body Design

Apex Volume




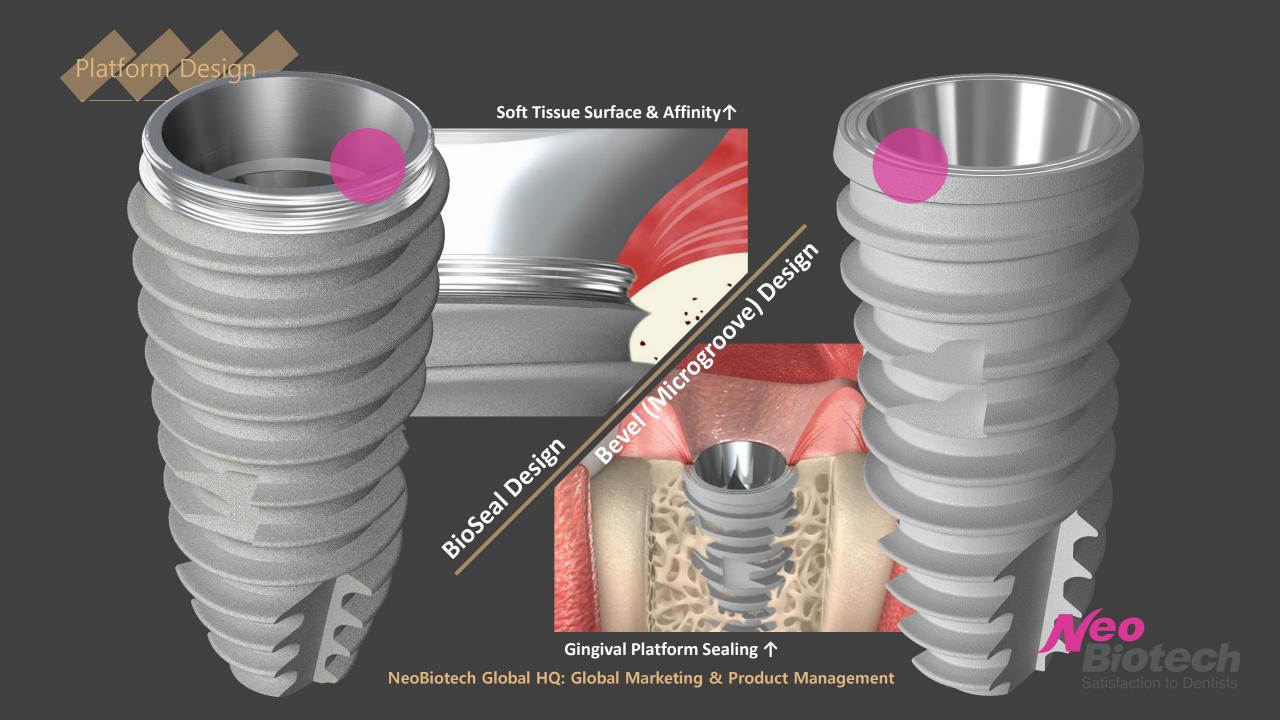

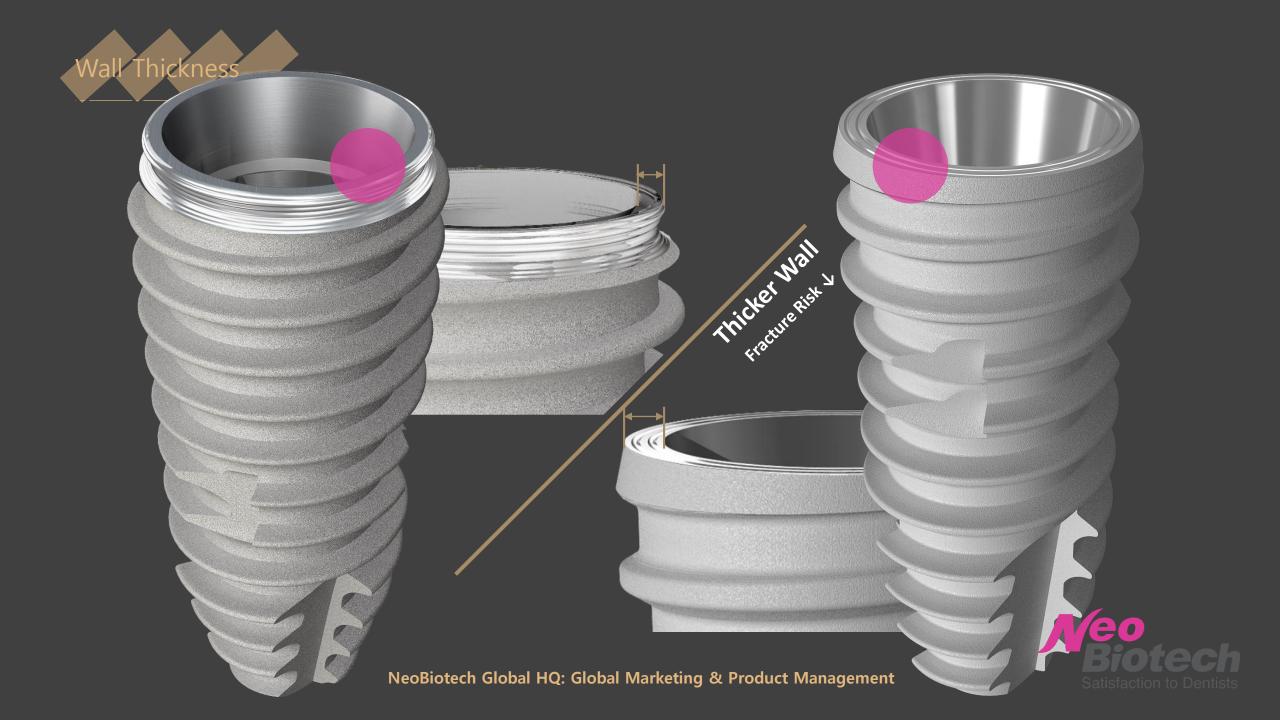



## Thicker wall of IS-III active than IS-III reduces fracture risk

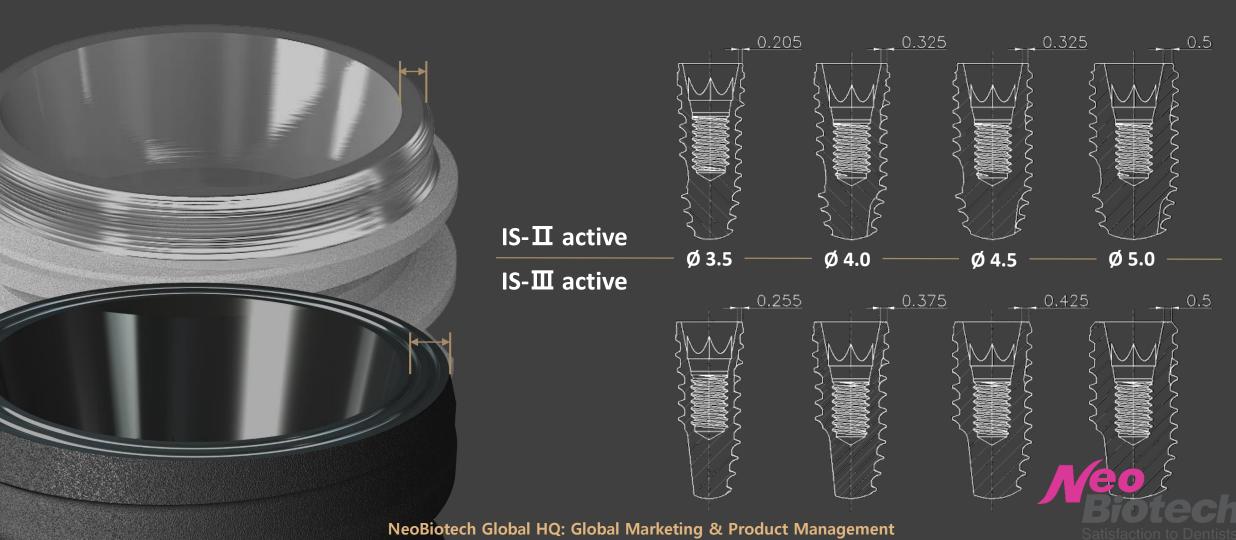

IS-**II** active


Self tapping ability is greater;

Fixture shows better performance with soft bone cases.

Overtorque might disrupt stability.


Thread on top 30% of fixture body is low;

For hard bone cases, tapping is required.

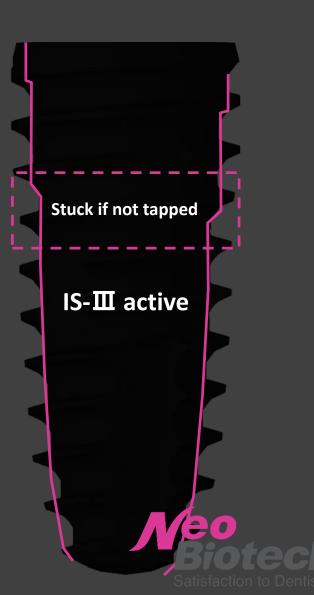




IS-**II** / IS-**III** Implantation

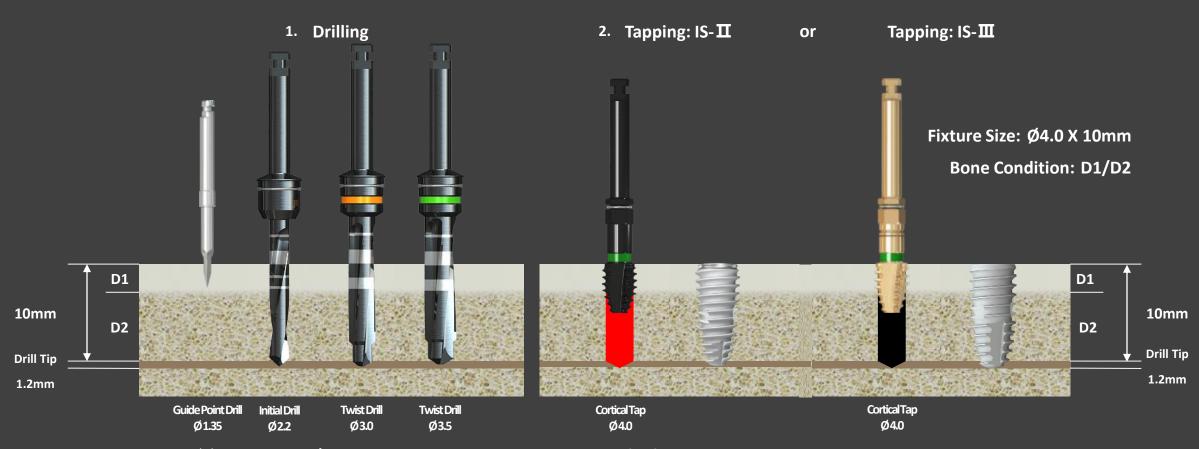

**X** Final drill is Ø3.0 surgical drill to implant in soft bone(D4) or to achieve more fixation at implant apex.



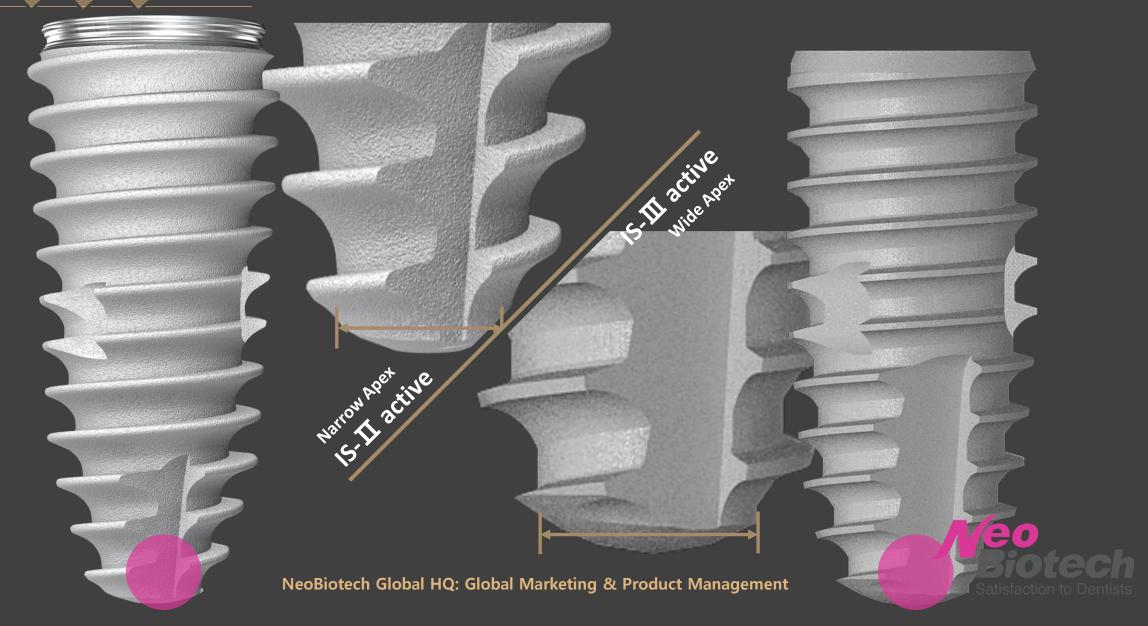

Apex: Volume



Apex of IS-**II** active is narrow, which shows excellent self-tapping ability.


For implantation case thru a narrow hole after sinus lift, the apex of IS-II active will give more stability.



Narrow Apex performs better penetration.

Bioteca Satisfaction to Pentists



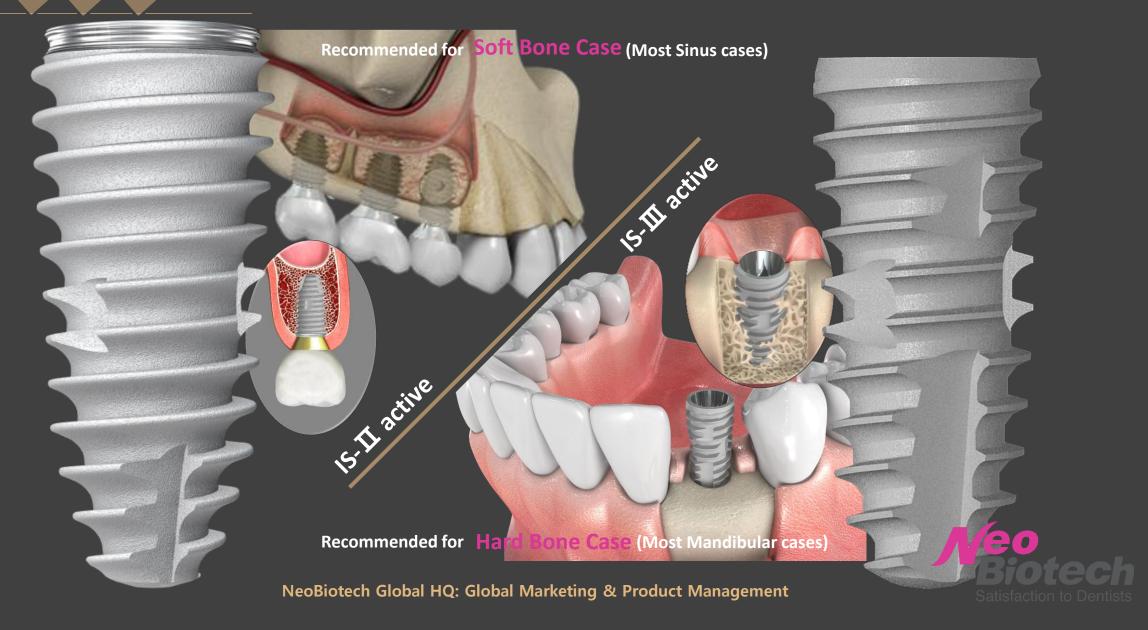

## Thanks

